

Results from Compound Optimization



Compound Optimization in MS and MS/MS

```

13:53:16: Optimizing S-Lens RF Amplitude for ion 336.18 m/z
13:53:22: Previous Setting = 150.00, New Setting = 101.22
13:53:22: Maximum Intensity = 1.68e+06
13:53:22: 479 % Improvement
13:53:23: Old Parent Mass: 336.176, New Parent Mass: 336.275
13:53:23: Optimizing collision energy at 1.5 mTorr
13:53:23: Waiting for the collision gas to stabilize
13:53:43: Finding the product ions of ion 336.3 m/z
13:54:22: Constructing the breakdown curve of ion 336.3 m/z
13:54:23: Product Ion: 247.15 Maximum Intensity: 1.49e+06
13:54:25: Product Ion: 179.06 Maximum Intensity: 8.72e+05
13:54:27: Product Ion: 191.05 Maximum Intensity: 7.93e+05
13:54:29: Product Ion: 158.17 Maximum Intensity: 4.37e+05
13:54:30: Product Ion: 205.09 Maximum Intensity: 1.36e+05
13:54:32: Product Ion: 102.11 Maximum Intensity: 4.58e+05
13:54:34: Product Ion: 144.10 Maximum Intensity: 2.58e+05
13:54:36: Product Ion: 156.10 Maximum Intensity: 3.55e+05
13:54:36: Collision Energy Optimization Results:
13:54:36: Product Ions (m/z) Coll. Energy (v)
13:54:36: 247.15 27
13:54:36: 179.06 43
13:54:36: 191.05 28
13:54:36: 158.17 26
13:54:36: 205.09 39
13:54:36: 102.11 39
13:54:36: 144.10 56
13:54:36: 156.10 43
13:54:38: Finish compound optimization
    
```

Comments:


Signature: _____

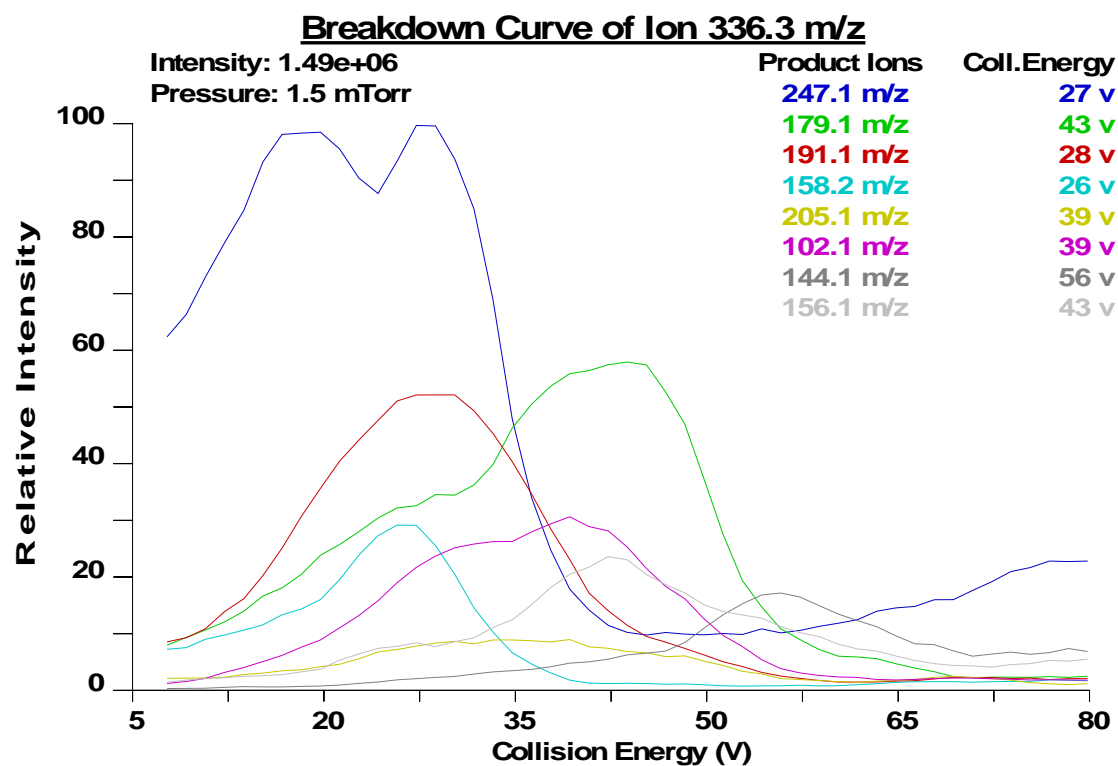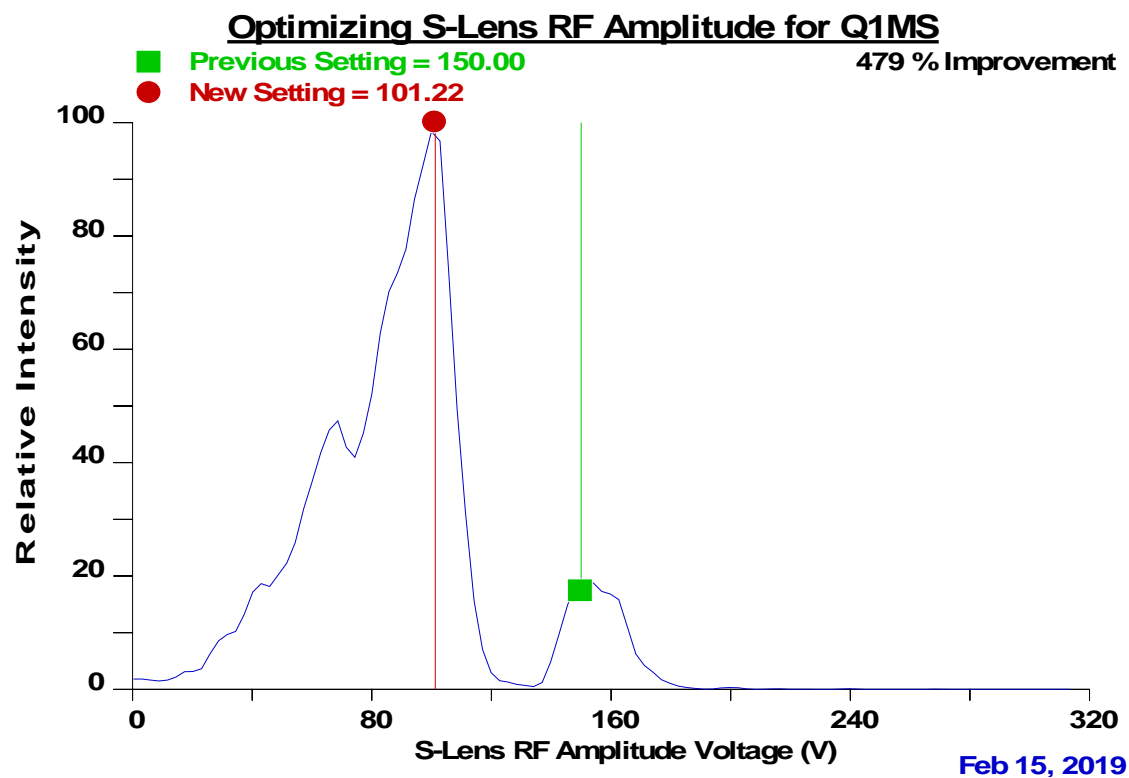


Signature: _____

TSQ Tune - Tune and Calibration

SRIGAMP DAC Scan #A: 19 Peak 99.64; 1.76e+006

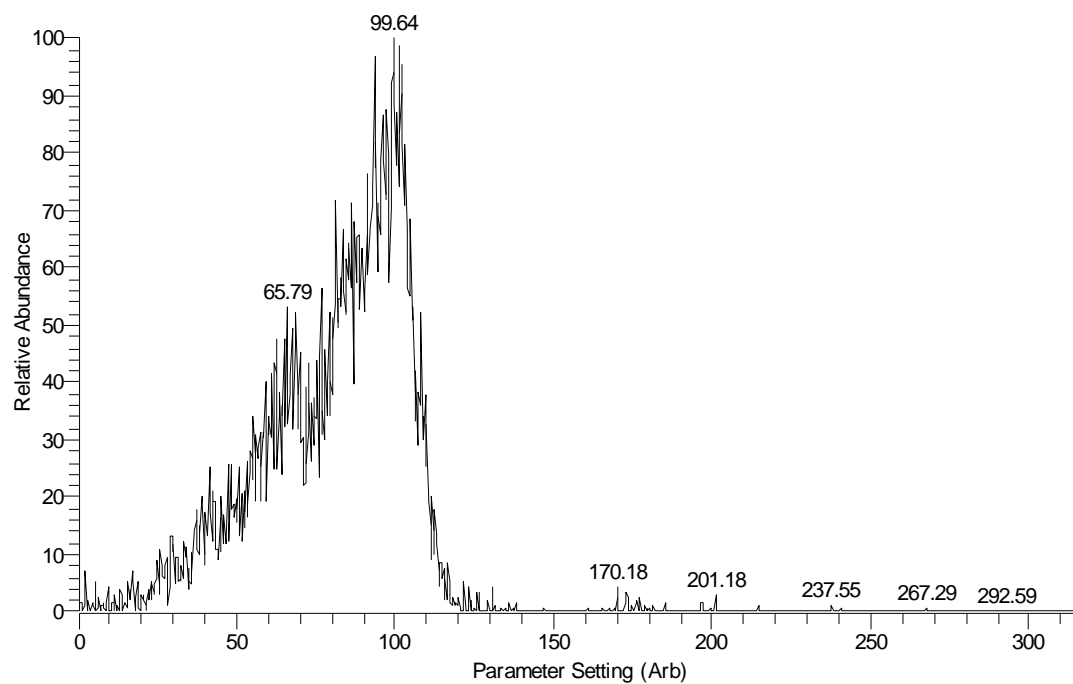

S#: 21270 FULL: PRO: 336 CT: 0.51 #A: 9 5.34e4

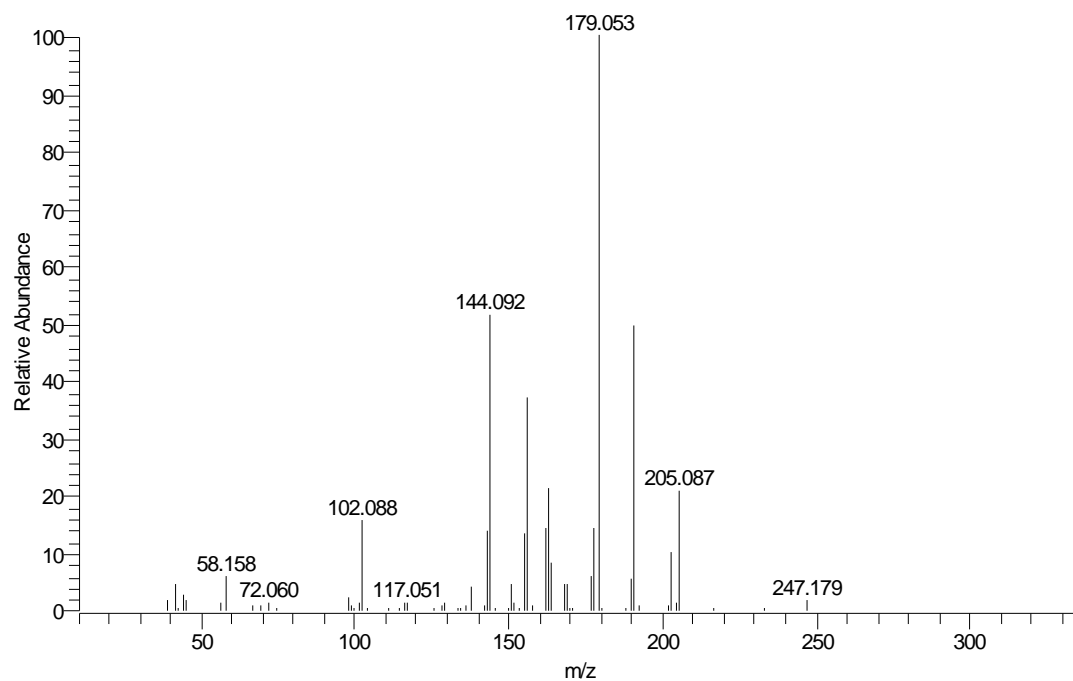

Signature: _____