

Results from Compound Optimization



Compound Optimization in MS and MS/MS

```

13:49:45: Optimizing S-Lens RF Amplitude for ion 336.18 m/z
13:49:50: Previous Setting = 200.00, New Setting = 105.33
13:49:50: Maximum Intensity = 1.41e+06
13:49:50: Significant Improvement
13:49:51: Old Parent Mass: 336.176, New Parent Mass: 336.332
13:49:51: Optimizing collision energy at 1.5 mTorr
13:49:51: Waiting for the collision gas to stabilize
13:50:12: Finding the product ions of ion 336.3 m/z
13:50:49: Constructing the breakdown curve of ion 336.3 m/z
13:50:50: Product Ion: 247.12 Maximum Intensity: 7.78e+05
13:50:52: Product Ion: 158.16 Maximum Intensity: 1.48e+05
13:50:54: Product Ion: 102.10 Maximum Intensity: 8.30e+04
13:50:56: Product Ion: 179.06 Maximum Intensity: 2.16e+05
13:50:58: Product Ion: 191.07 Maximum Intensity: 1.22e+05
13:50:59: Product Ion: 69.19 Maximum Intensity: 8.71e+03
13:51:01: Product Ion: 58.09 Maximum Intensity: 1.92e+04
13:51:03: Product Ion: 72.06 Maximum Intensity: 9.03e+03
13:51:03: Collision Energy Optimization Results:
13:51:03: Product Ions (m/z) Coll. Energy (v)
13:51:03: 247.12 20
13:51:03: 158.16 20
13:51:03: 102.10 27
13:51:03: 179.06 37
13:51:03: 191.07 31
13:51:03: 69.19 31
13:51:03: 58.09 42
13:51:03: 72.06 51
13:51:05: Finish compound optimization
    
```

Comments:


Signature: _____

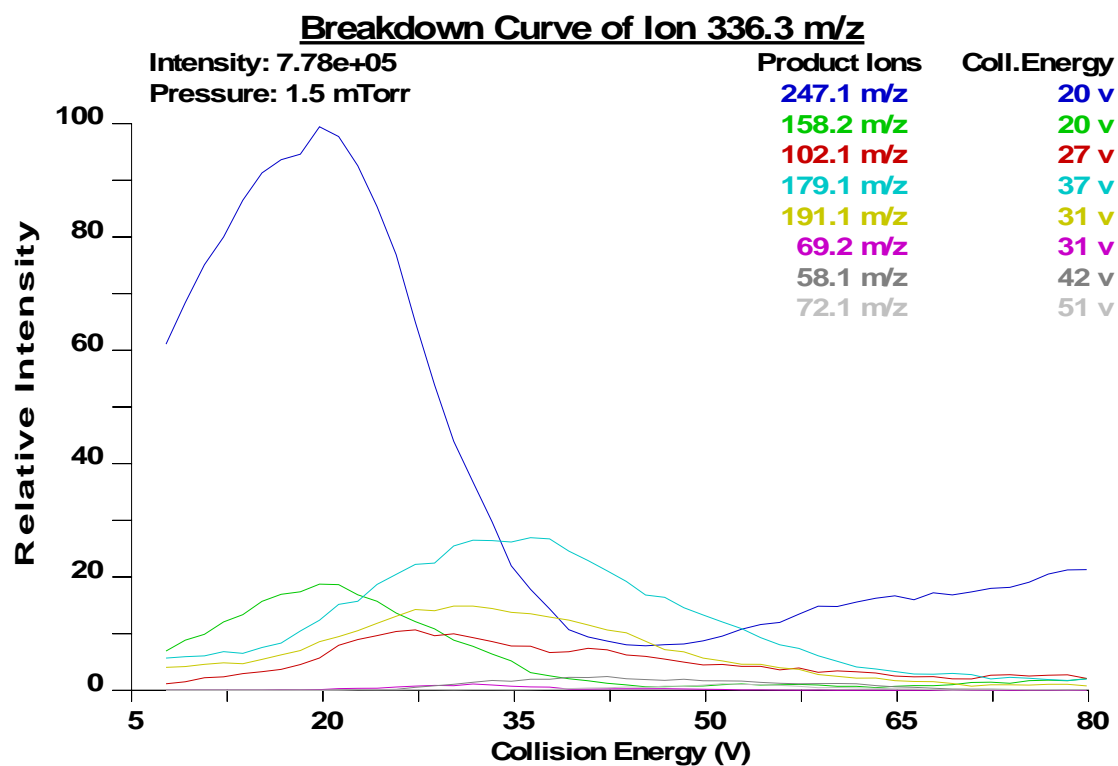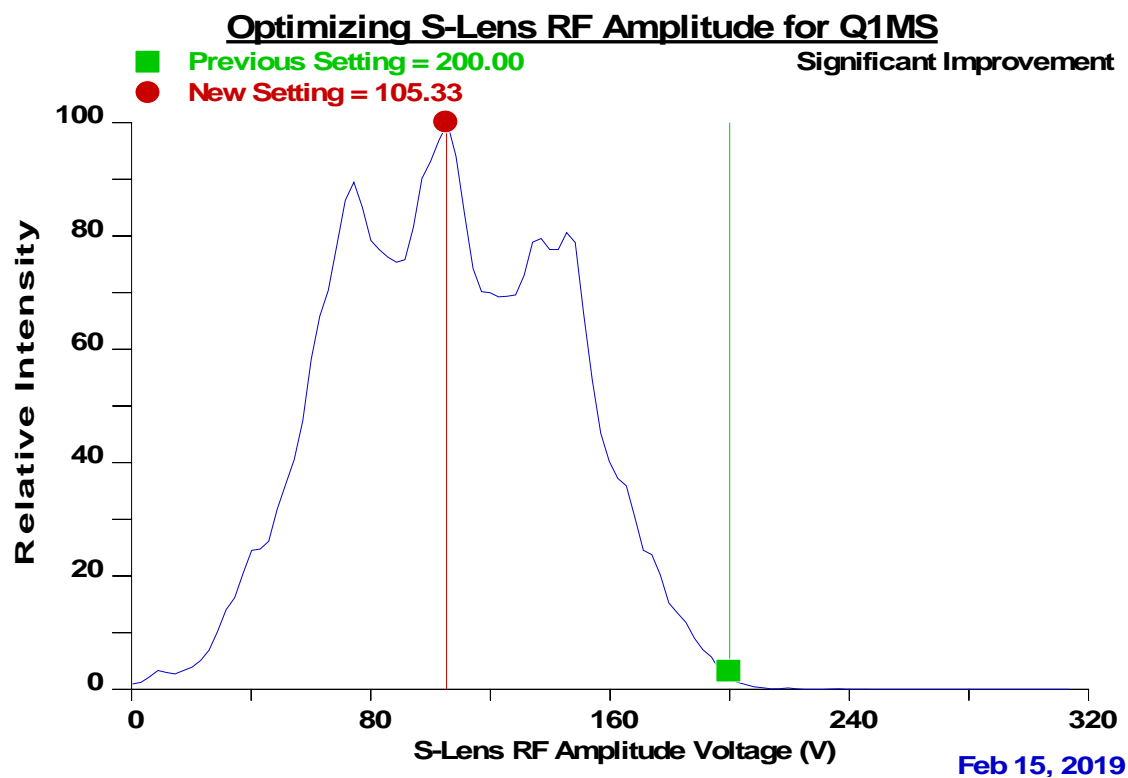


Signature: _____

TSQ Tune - Tune and Calibration

SRIGAMP DAC Scan #A: 19 Peak 103.75; 1.48e+006

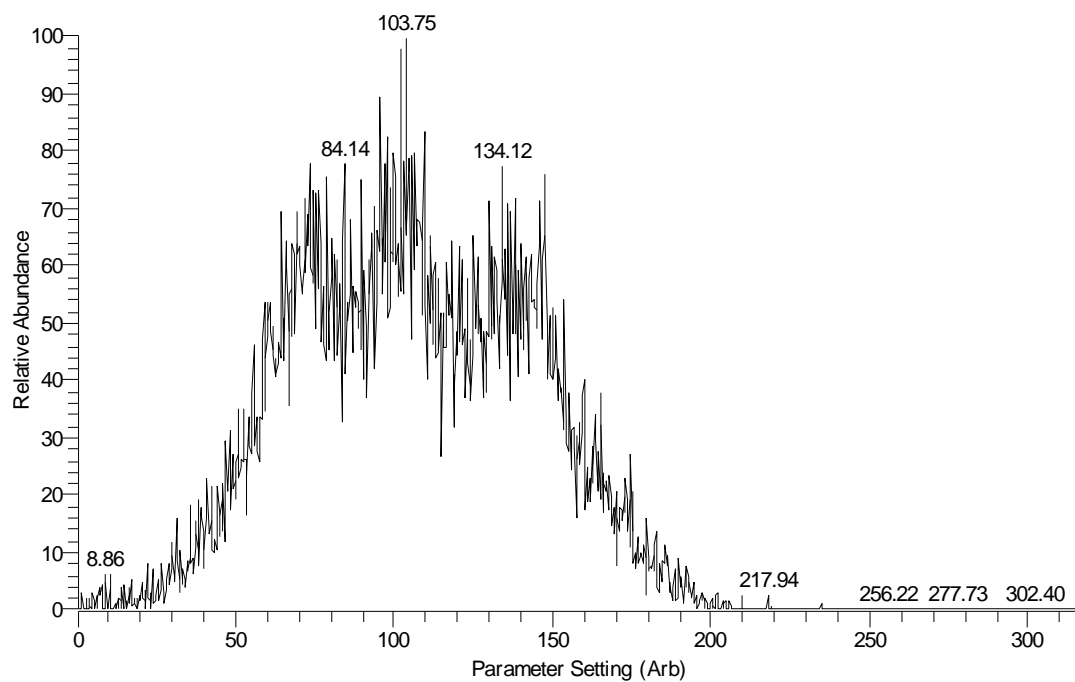

S#: 20718 FULL: PRO: 336 CT: 0.51 #A: 9 7.11e3

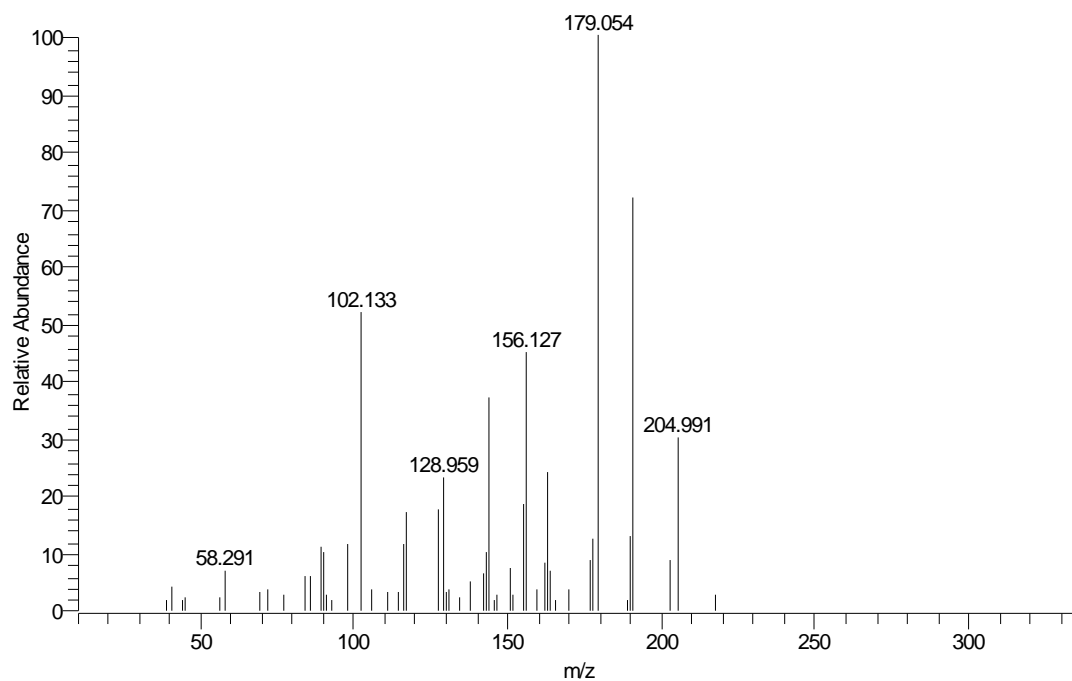

Signature: _____